


## PA-E2000-4

Polarization Analyzer adapters for fiber cables without optics

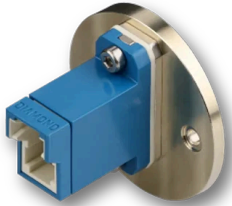

### FEATURES

Polarization Analyzer adapters for fiber cables without optics.

- E2000 APC
- Inclined coupling axis

This adapter can also be used with fiber connectors of type E2000 PS. Please note, that in this case the iris diaphragm of the Polarization Analyzer cannot be used to reduce laser power as the outgoing beam is too small. Please keep this in mind when analyzing radiation with high power.

## FAQ

### Polarization Analyzer accessories

**There are two adapter options - one without integrated optics and one with integrated optics. What is the best option for me?**

For standard applications we recommend to use an adapter without optics since any additional optics may show birefringence and so may falsify the measurement result.

This is a printout of the page <https://sukhamburg.com/products/details/PA-E2000-4-0> from 5/7/2024

## CONTACT

For more information please contact:

**Schäfter + Kirchhoff GmbH**

**Kieler Str. 212**

**22525 Hamburg**

**Germany**

**Tel: +49 40 85 39 97-0**

**Fax: +49 40 85 39 97-79**

**[info@sukhamburg.de](mailto:info@sukhamburg.de)**

**[www.sukhamburg.com](http://www.sukhamburg.com)**

## LEGAL NOTICE

**Copyright 2020 Schäfter+Kirchhoff GmbH. All rights reserved.**

Text, image, graphic, sound, video and animation files and their arrangement on Schäfter+Kirchhoff GmbH webpages are protected by copyright and other protective laws. The content may not be copied for commercial use or reproduced, modified or used on other websites. [\[more\]](#)